

ThermaKnife™

Hot Tipped Tray Trimming Knife



An effective and time-saving knife for fast and easy trimming of bleaching trays, night guards, mouthguards and other thermoplastic materials.

Featuring an A and B switch, the ThermaKnife has the ability to support two different pens at one time or support a pen and our newly available ThermaPlier. This allows the user to quickly switch back and forth on a trimming job for efficiency and accuracy. Simply choose the pen design that best suits your needs, set your temperature and complete your trimming job with ease. With an adjustable temperature control, you'll always be able to customize the process to your cutting speed and material thickness. Pens and plier can also be interchanged quickly and easily.




ThermaPlier NEW

Hot Tipped Thermal Pliers

Our new ThermaPlier hot tipped thermal plier attachment is great for reheating the critical areas of dental aligners and tooth positioners without causing adjacent plastic to distort. Fit and form are very important when dealing with aligners, as they function as a guide to proper tooth position, so trust in a tool to get the job done right! It can also be used in orthodontic appliances and mouthguards.



*Optional Accessory

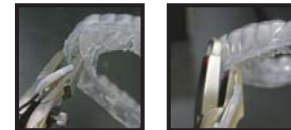

Features/Benefits:

- Fully adjustable temperature control for precise cutting action
- Cutting tip heats ultra-fast to save time and eliminate waiting
- Lightweight pen features an ergonomic design and sharp, permanent tip for improved comfort & control
- A & B Selection switch allows for the use of two different Pens or a Pen and ThermaPlier to reduce time and increase trimming productivity
- Fast and easy pen replacement (via snap connector at back of pen) for added convenience
- Convenient pen rest provides added safety when unit not in use
- Small portable size: Console 5.75" W x 5" D x 4" H (incl. pen rest) Pen: 6" L x .75" W (at foam grip), 2.5 lb total console/pen weight

Ordering Information:

Item No.	Description		
80500	ThermaKnife Intro Kit, 120 V AC (Includes Knife-Edge Pen)	80530	Replacement Pen with Waxing Tip
80510	Replacement Pen with Concave Tip	80610P	ThermaPlier
80520	Replacement Pen with Knife-Edge Tip		

Contact your dealer today, or jeannieb@mailbdm.com and we will be happy to assist you! ThermaKnife and ThermaPlier are Made in the USA!

Buffalo Dental Manufacturing Co., Inc.

159 Lafayette Dr., PO Box 678, Syosset, NY 11791-0678 | Tel: 516.496.7200 or 800.828.0203
 Fax: 516.496.7751 | Email: sales@mailbdm.com | www.buffalodental.com